KÖNIG
ELECTRONIC

HC-PL10-20-30

MANUAL (p. 2)
MASSAGE PILLOW

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Română

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ENGLISH


This product can effectively massage muscles, release arteries and veins and increase the blood circulation. It relieves stiffness, tiredness and paralysis while it can be used wherever: at home or at the office.

Product Features

1. 6 massage modes (Ref. to table 1)
2. Automatically switches off after 20 minutes
3. Improve the blood circulation of your head and neck for a better oxygen supply to the brain
4. The massaging device can be installed inside a pillow which can adapt to the different postures of the users and provide good support. This supportive function minimizes the harm caused by pressure and protects the backbone.
5. This pillow is designed to fit a human's head curve ergonomically.

Table 1

Mode	Massage Type
MODE 1	Continuous gentle massage
MODE 2	Continuous strong massage
MODE 3	Changes from strong massage to stop massaging (Changes 70 times within 1 minute)
MODE 4	Alternates gentle and strong massages. Changes between the two modes 60 times in 1 minute.
MODE 5	Circulate massage mode 1,2,3 and 4 for around 3 seconds each
MODE 6	Circulate massage mode 1,2,3,4,and 5 for around 45 seconds each

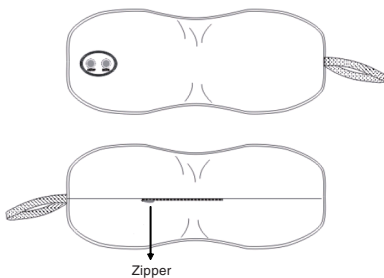1 Operation

- 1.1 Press the ON/OFF button to switch on the massage pillow and activate the massage mode.
- 1.2 Press the MODE button to select the different massaging modes.

2 Battery Replacement

When the massage pillow's vibration is slower or stops completely, the battery is low. Please replace the battery with 2pcs of AA batteries.

- 2.1 Unzip the back of the pillow and take out the battery compartment.
- 2.2 Press on the battery cover slightly and slide it down according to the arrow's direction on the cover.
- 2.3 Replace the old batteries by the new ones.
- 2.4 Close the battery cover. Put it back inside the massage pillow and zip the massage pillow again.

3 Attention

- 3.1 Please put the massage pillow away from children.
- 3.2 Do not use or store the massage pillow near magnetic fields and keep it away from high temperatures and moisture.
- 3.3 Please take out the batteries and store the massage pillow properly when it is not going to be used for a long period.
- 3.4 People without the ability to express themselves, such as babies, are forbidden to use this product.
- 3.5 Do not throw the batteries into a fire.

Declaration of conformity

We,
Nedis B.V.
De Tweeling 28
5215MC 's-Hertogenbosch
The Netherlands
Tel.: 0031 73 599 1055
Email: info@nedis.com

Declare that product:
Brand: König Electronic
Model: HC-PL10-20-30
Description: Massage pillow

Is in conformity with the following standards
EN55022:2006
EN61000-6-3:2001+A11:2004
EN61000-6-1:2001

And complies with the requirements of the European Union Directives.
2004/108/EC

's-Hertogenbosch, 2009, January 22

NEDIS B.V.
Postbus 2009 5201 CZ 's-Hertogenbosch NL
De Tweeling 28 5215MC 's-Hertogenbosch NL
Telefoon 073-599 96 41

Mrs. J. Gilad
Purchase Director

Safety precautions:

CAUTION
RISK OF ELECTRIC SHOCK
DO NOT OPEN

To reduce risk of electric shock, this product should **ONLY** be opened by an authorized technician when service is required. Disconnect the product from mains and other equipment if a problem should occur. Do not expose the product to water or moisture.

Maintenance:

Clean only with a dry cloth. Do not use cleaning solvents or abrasives.

Warranty:

No guarantee or liability can be accepted for any changes and modifications of the product or damage caused due to incorrect use of this product.

General:

Designs and specifications are subject to change without notice. All logos brands and product names are trademarks or registered trademarks of their respective holders and are hereby recognized as such.

Attention:

This product is marked with this symbol. It means that used electrical and electronic products should not be mixed with general household waste. There is a separate collections system for these products.

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Copyright ©
Copyright © (Autorská práva)